Allcem Core fulfilled three distinct functions (cementation of the post, core build up, and cementation of the crown) providing the work with greater agility.

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
INITIAL | INICIAL

CASE HIGHLIGHTS | DESTAQUES DEL CASO

After the removal of the crown and core, Top dam (gingival barrier) aids in the absolute isolation of the tooth with little coronal structure.

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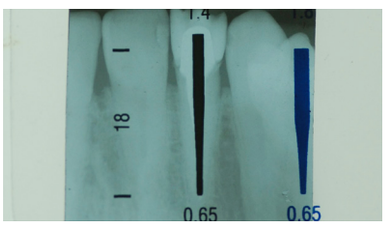
The choice of the post using the Whitepost template is very simple and is done by simply overlapping it on the periapical X-ray of the tooth.

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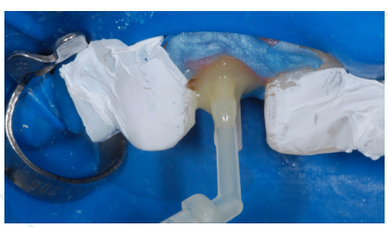Using the intracanal tip for self-mixing of the Allcem Core resin cement, it is possible to inject cement inside the already-prepared canal simply and quickly. After the installation of the fiberglass post, the same cement can be used to build the core.

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Core made with Allcem Core, ready to receive the crown.

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EN

Reconstruction of endodontically treated teeth with **fiberglass post and lithium disilicate crown**

AUTHORS: Arcelino Farias Neto, Monalisa M. Queiroz Jales, Simone L. Pinheiro de Sousa and Cláudia T. Machado Cunha

Male patient, 57.

COMPLAINT

Crown of the upper right lateral incisor presenting a color different from the other teeth.

INITIAL APPEARANCE

Metal ceramic crown with whitish and opaque appearance, with unsatisfactory fit.

CLINICAL/X-RAY EXAM

Tooth 12 presented a cast metal core with an intraradicular length less than the ideal, and a metal ceramic crown of a color incompatible with the other teeth and deficient marginal fit. The endodontic condition of the tooth was satisfactory on the periapical X-ray and the patient presented no negative clinical symptoms.

TREATMENT

It was proposed to remove the crown and the metal core. A fiberglass post was installed in the root to support the core, over which the lithium disilicate crown would be seated.

> "Allcem Core cement was used in the three cementation steps, which allowed the work to be done faster."
>
> *DR. ARCELINO FARIAS NETO*

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FGM MATERIALS USED

DUAL-CURING RESIN CEMENT
Allcem Core (color A2)

ADHESIVE FOR ENAMEL AND DENTIN
Ambar

FIBERGLASS POST
Whitepost (DC 0.5)

GINGIVAL RETRACTION CORD
Pro-Retract (number 00)

DISPOSABLE MICRO APPLICATOR
Cavibrush long

2% CHLORHEXIDINE SOLUTION
Clorhexidina S

37% PHOSPHORIC ACID ETCHANT
Condac 37

10% HYDROFLUORIC ACID ETCHANT
Condac Porcelaina

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