
X-ray appearance of tooth 46. Observe the correct fit of the fiberglass post in the root canal and fill core made, giving the tooth a final sealing and fracture resistance, ready to receive the prosthetic work.

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
INITIAL | INICIAL

CASE HIGHLIGHTS | DESTAQUES DEL CASO

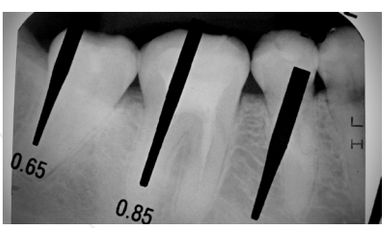In the initial planning, the template of the Whitepost system allows for predicting the best post indicated for the case.

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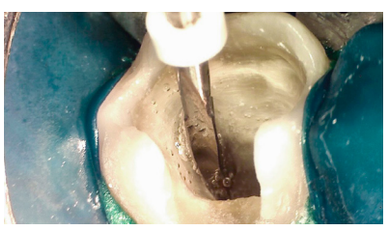The Whitepost DC1 drill used when instrumenting the root canal confers the ideal tapered shape to the canal for the post to fit.

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Intraoperative X-ray demonstrating the correct procedure for deobturation of the distal canal.

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
Cementation of the fiberglass post (Whitepost DC1, FGM). Allcem Core dual-curing resin cement (FGM) offers a tip that allows the cement to be applied inside the prepared canal, which speeds up the procedure. This way, the same cement allows for immediate making of the core of the cavity.

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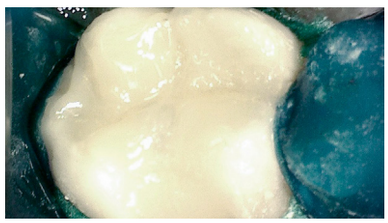Final clinical appearance of tooth 46 restored with Opallis composite resin (FGM).

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EN

Socket–shielding technique: installation of **direct posts** as retainers **for extensive restoration**

AUTHORS: J. Edgar Valdivia C., Hector Caballero Flores and Manoel E. de Lima Machado

Male patient, 27.

COMPLAINT

Lower molar with pain when chewing and with extensive restoration.

INITIAL APPEARANCE

Initial X-ray appearance of tooth 46 showing the need for endodontic re-treatment.

CLINICAL/X-RAY EXAM

Tooth 46 showed extensive occusal/vestibular restoration requiring replacement, and the X-ray showed the need for endodontic re-treatment and intraradicular reinforcement.

TREATMENT

Endodontic re-treatment was planned, aiming for installation of a fiberglass post and restoration with composite resin in a single session. For this purpose, at the time of deobturation and preparation of the distal canal, the Whitepost DC1 drill (FGM) was used, leaving the canal prepared for the installation of the Whitepost DC 1 fiberglass post (FGM). It was immediately proceeded with cementation of the post and core of the cavity (Allcem Core, FGM) and completion of the restoration with composite resin (Opallis, FGM).

> "The combination of Whitepost fiberglass posts and Allcem Core resin cement allowed for the simultaneous radicular and coronal restoration to be done in a practical and safe manner."
>
> *DR. J. EDGAR VALDIVIA CARDENAS*

FGM MATERIALS USED

DUAL CURING RESIN CEMENT
Allcem Core

FIBERGLASS POST
Whitepost (DC 1)

NANOHYBRID COMPOSITE RESIN
Opallis

37% PHOSPHORIC ACID ETCHANT
Condac 37

GINGIVAL BARRIER
Top dam

SILANE
Prosil

DISPOSABLE MICRO APPLICATOR
Cavibrush regular

ADHESIVE FOR ENAMEL AND DENTIN
Ambar

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