FGM

Allcem CORE

2797

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EN Instructions Manual

Dual-curing resin cement for core build-up and cementation of intraradicular posts and crowns.

Professional use only

Read all the information in this document carefully before using the product. Keep the document for later reference at least until the product is finished or until there is no more interaction between the product and its last patient.

Product description

Allcem Core is a dual-curing resin cement specially conceived for the adhesive cementation of intraradicular posts, crowns and core buildup (vital and non-vital teeth). The product is commercialized in a double-body syringe, which facilitates obtaining the correct mixing proportions of the two pastes (1:1), and eases clinical procedures. By using a self-mixing tip with an applying nozzle (intracanal tip), it is possible to apply the cement inside the root canal and easily cement the post. The same tip can be used to add material to the crown region of the post to build up the core. Allcem Core is a composite made off methacrylate monomers such as BisGMA, BisEMA and TEGDMA, inorganic load, photoinitiators, coinitiators, catalysts and pigments. The cement associates the necessary flow for applications in canals with the ideal thixotropy for core build-ups. The cement has mechanic strength comparable to that of resin composites, allowing for building up the core with safety.

Presentation forms

· Dual-body syringe (base + catalyst) with NET WT. 0.21 Oz (6g) and 8 self-mixing tips with applicator nozzle. Available shades: A1, A2, A3 (Universal) and Opaque Pearl.

· Packaging with 20 disposable self-mixing tips with applicator nozzle.

Basic Composition

Base paste: methacrylate monomers such as BisGMA, BisEMA and TEGDMA, camphorquinone, coinitiators, barium-aluminum-silicate glass microparticles and particles of silicon dioxide, inorganic pigments and preservatives. Catalyst paste: methacrylate monomers, dibenzoyl peroxide, stabilizers and barium-aluminum-silicate glass. The mixture of base paste and catalyst

pastes at proportion (1:1), presents approximately 82% load in weight.

Product indication

Allcem Core is a dual-curing radiopaque permanent resin cement system indicated for the cementation of indirect restorations to the tooth structure. It can be used for permanent cementation of;

· Porcelain fused-to-metal crowns and bridges;
· Metal crowns, bridges, inlays and onlays (high noble, noble, and base metals);
· Crowns and bridges with minimal tooth structure;
· Maryland bridges (resin bonded bridges);
· All ceramic/porcelain and pre-cured composite crowns, bridges, inlays and onlays;
· Endodontic posts;
· Adhesive core build-up of vital and non-vital teeth

Note: For cementation of laboratory-fabricated veneers, use Allcem Veneer APS (FGM).

Warnings and Precautions

During product handling

· Allcem Core cement contains acrylate resins. Avoid use of this product on patients with known acrylate allergies. To reduce the risk of allergic response, minimize exposure to these materials. In particular, avoid exposure to uncured resins;
· Use of protective gloves and a no-touch technique is recommended;
· If skin contact occurs, wash skin with soap and water;
· Acrylates may penetrate commonly-used gloves. If cement contacts glove, remove and discard glove, wash hands immediately with soap and water and then re-glove;
· If accidental contact with eyes or prolonged contact with oral soft tissue occurs, flush with large amounts of water. If irritation persists, consult your physician.

During product application

· If there is pulp exposure, protect the area with a base/liner material;
· Do not use eugenol-based provisional cements before using Allcem Core;
· Allcem Core is not compatible with single-step self-etching adhesives, except Ambar Universal APS (FGM), which can be used according to its protocol;
· Light cure the resin cement for at least 40 seconds on each margin, with light curing unit (halogen light or LED) with power of 450mW/cm² or more and blue light emission spectrum of (400-500 nm). Wear protective goggles when operating a light-curing unit;
· Discard the self-mixing tip after use;
· Be sure to place back the cap after use.

Contraindication


The product should not be used in patients allergic to any of the substances contained in the formula.

Side effects

Systemic side effects are not known to date. In allergic individuals, allergic reactions may occur.

Instructions for Use

Table 1 : Coupling the self-mixing tip to the double-body syringe.

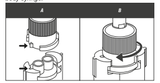a) After removing the cap of the double-body syringe, align the triangular-shaped protruding mark at the base of the self-mixing tip with the mark at the top of the syringe. The alignment and coupling of the two parts must occur in this position

b) After coupling, spin (clockwise) the self-mixing tip ¼ of a turn for the correct final coupling. Make sure the coupling is perfect and ready for use.

Obs.: The self-mixing tip with applicator nozzle is disposable (a) After removing the cap of the double-body syringe, align the triangular-shaped protruding mark at the base of the self-mixing tip with the mark at the top of the syringe. The alignment and coupling of the two parts must occur in this position

b) After coupling, spin (clockwise) the self-mixing tip ¼ of a turn for the correct final coupling. Make sure the coupling is perfect and ready for use.

Obs.: The self-mixing tip with applicator nozzle is disposable (use only once).

Tooth preparation

1. Isolate the tooth and use a retraction system when the margin of the preparation is sub-gingival to make sure no contamination occurs in the preparation.

2. For cementing intraradicular posts, the canal must be prepared following the instructions of the system used (Whitepost, FGM). Once prepared, the canal must be etched with phosphoric acid at 37% during 15 seconds then washed profusely with water and dried, but left slightly moist for the adhesive application. In the case of core build-up (vital or non-vital teeth), the acid etching is the same on the dental remains preparing it to receive the adhesive.

3. Apply the adhesive following the manufacturer´s instructions and light cure it. It is recommended the use of two or three-step adhesives. Dual-curing resin cements are incompatible with self-etching adhesives (single-step), except Ambar Universal APS (FGM), which can be used according to its protocol.

Cement characteristics

Table 2: Times.

a) Mixing time spatula (not applicable when using auto mixing tips).		*20-30"*
b) Working time. Application on the canal or increments.		*1'30"*
c) Chemical curing time.		*6-7'*
d) Photo-activation. Per increment.		*40"*

Cementing intraradicular posts and building up the core.

1. After an adequate preparation of the tooth (isolation, partial removal of gutta-percha, acid etching and careful adhesive application), post the (Whitepost, FGM) must be prepared as needed. A structure cut is indicated by using drills and silanization of its surface with silane according to the instructions of the manufacturer.

2. For the insertion of Allcem Core in the canal, there are two suggested procedures:

a) Manual: mix the base paste and the catalyst paste in proportion (1:1) on a glass surface or a mixing pad for 20-30 seconds. After obtaining a homogeneous mass, the mixture must be applied in the canal with fine instruments (e.g.: straight probe) or lentulo. The cement can also be applied to the post surface and then put in place.

b) Self-mix: Allcem Core has a self-mixing tip with applicator nozzle, adequate for the direct insertion of the cement to the interior of the canal. In this case, make sure the self-mixing tip and the double-body syringe are correctly coupled.

3. Position the post in the root canal. For building up the core, there are three suggested alternatives:

a) Use excess cement that overflows from post cementation and, by using spatulas, form the core. Light-cure it after the cementation.

b) The crown part of the post can receive cement increments by using the same syringe used in the intracanal insertion. These increments can be light cured until the completion of the final shape of the core.

c) As an option, the professional can use pre-manufactured acetate matrices with the form of the core. In this case, it is recommended to fill out the mold with the cement and adapt it to the crown part of the post and then light-cure it independently of the technique used, the increments must be light-cured for 40 seconds each.

4. After adequate light curing of the core, proceed with the final forming of it using disks, sanding, drills.

Building up cores in vital or non-vital teeth without an intraradicular post:

1. To build up the core in vital or non-vital teeth which did not receive intraradicular posts, perform acid etching and adhesive application following the adhesive technique previously described.

2. The core build-up can be executed in successive increments of cement by using either the self-mixing tip or a spatula for composites.

3. Each increment must be light cured for 40 seconds (*).

4. After adequate light curing of the core, proceed with the final forming of it using disks, sanding, drills.

*Considering a light-curing unit of minimum power of 450mW/cm² and spectrum of blue light (400 a 500 nm).

Cementation of crowns of pure ceramic, metal-ceramic or laboratory composites:

1. Prosthetic pieces must receive previous treatment depending on the material.

a) Metallic or metal-ceramic: it is recommended the jet sanding with abrasive particles (in the laboratory) to create the necessary surface roughness to improve the anchoring of the piece. It can also be considered the use of a primer for metal (following the manufacturer´s instructions) to improve the retention strength.

b) Pure ceramic (etchable) or laboratory composites: it is suggested the acid etching of the internal surface with hydrofluoridric acid at 5% and 10% and posterior application of silane following the corresponding product instructions.

2. After the treatment of the piece, prepare the tooth following the adhesive protocol previously described.

3. Select the adequate shade of Allcem Core for the cementation of the piece. To handle the product, it can be dispensed on a mixing pad (base and catalyst) in equal proportions and proceed to its homogenization for 20-30 seconds. Alternatively, there is the option of using the self-mixing tip with applicator nozzle to dispense the product directly on the piece.

4. Apply a fine layer of cement on the internal part of the piece and carefully position it in place confirming an adequate overflow of cement over the margins of the preparation.

5. Keep the piece in position and start the cleaning of the excess cement after 2 minutes of the initial mixture of the pastes. To avoid the effect of inhibited layer, a layer of solid vaseline can be applied on the cement excesses. In case the excess of cement was removed immediately after the positioning of the piece, the margins must be light cured to avoid the inhibited layer effect due to oxygen.

6. After removing the excess of cement, light cure every surface and all margins for 40 seconds (*) or allow the chemical curing to continue for 7 minutes counting from the initial mixing of the pastes and considering also the room temperature.

Note 1: For porcelain and resin composite pieces, the margins can be light cured for 40 seconds (*) to ensure the immediate fixation of the piece.

Note 2: After light curing the restorations and the margins, remove the excess cement by using a scalpel. The use of rotary instruments must be avoided before 10 minutes from the initial mixing of the pastes, depending also on room matériaux.

7. Instruct the patient not to exert force on the restorations for at least 15 minutes after cementation.

8. Proceed with the polishing and finishing procedures.

* Considering a light-curing unit with a minimum power of 450mW/cm² and spectrum of blue light (400 a 500 nm).

Conservation and storage

Syringe must be immediately closed after each use.

Do not expose the product to high temperatures or intense light.

Store the product at temperatures between 5°C and 25°C / 41°F and 77°F.

Expiration date is printed on the packaging.

Warnings

All dental products should be used only by dentists, dental prosthetics and authorized laboratories. Do not use the product after it is expired. When disposing the product, follow the legislation of your country. Keep out of the reach of children.

This material was manufactured for dental use only and must be manipulated according to the instructions for use. The manufacturer is not responsible for any damage caused by misuse or improper manipulation of the product. The user is obligated to prove (before using the product and under his/her responsibility) that this product is compatible with the desired use, especially when such use is not indicated in the instructions for use. Data description does not constitute any type of warranty and therefore is not binding.

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