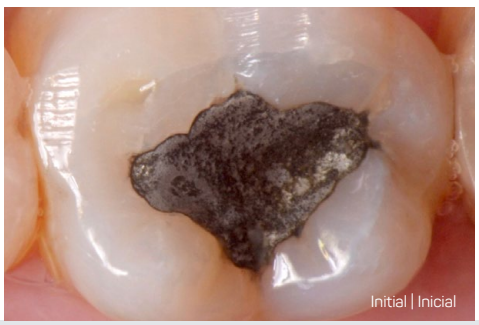
EN One-step direct restoration

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40-year-old female patient.

COMPLAINT

Sensitivity to thermal stimulation in a tooth with an old restoration.

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EN CLINICAL/RADIOGRAPHIC EXAM

Presence of amalgam restoration in the occlusal face of tooth 46 with an oxidized porous surface and absence of anatomy.

TREATMENT

Removal of the old restoration and preparation of a new restoration with bulk-fill composite (Opus Bulk Fill APS) in a single increment. The chosen material provides greater agility associated with aesthetics.

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

2. Isolated tooth.
2. Diente aislado.

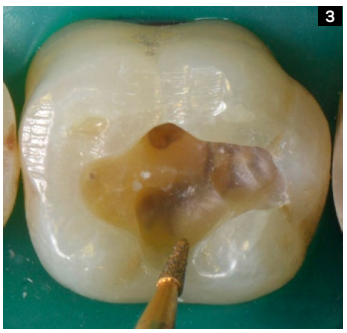
3. Cavosurface angle smoothing.
3. Alisado del ángulo cavo superficial.

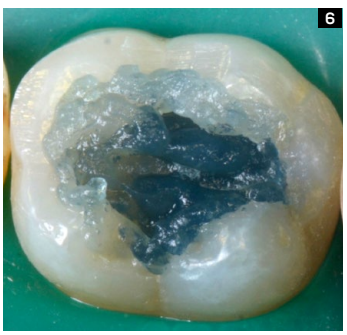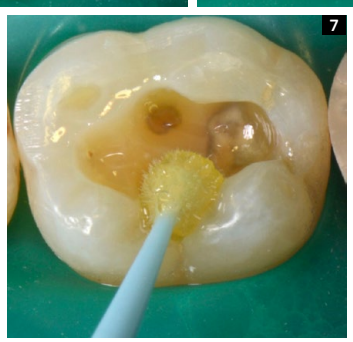4. Cavity ready.
5. Enamel acid etching for 30 seconds (Condac 37).
6. Dentin acid etching for 15 seconds (Condac 37).
7. Application of Ambar adhesive.

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8. Light-cured adhesive.
8. Adhesivo fotocurado.

9. Insertion of Opus Bulk Fill APS composite in a single increment.

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Agility when performing posterior restorations in large increments is a feature that brings many benefits to the dentist's daily routine. Opus Bulk Fill APS is my composite of choice, because it has exceptional consistency and aesthetics, in addition to the APS system that provides greater safety regarding the cure depth.

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Profa. Dra. Dayse Lúcia Otero Amaral

10. Restoration finished with isolation

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11. Restoration finished without isolation after polishing with Diamond felt disks and Diamond R polishing paste.

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LIGHT-CURING BULK-FILL COMPOSITE
Opus Bulk Fill APS

RESINA FOTOCURABLE TIPO BULK FILL
Opus Bulk Fill APS

5TH GENERATION LIGHT-CURING ADHESIVE
Ambar

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