EN Direct restoration of posterior teeth in a practical and effective way

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29-years old female patient.

COMPLAINT

Pain stimulated by cold stimulus in the left inferior region.

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EN CLINICAL/RADIOGRAPHIC EXAM

Teeth number 36 had a restoration with marginal infiltration besides deficiency in shade and form.

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TREATMENT

The restoration was redone with the easy-to-apply Opus Bulk Fill composite, which has low shrinkage stress. Tooth number 36's restoration was performed in two steps with the same composite. Initially, the proximal cavity was closed and after that, the occlusal cavity was filled. Tooth number 37 had a class I cavity which was restored in a single step.

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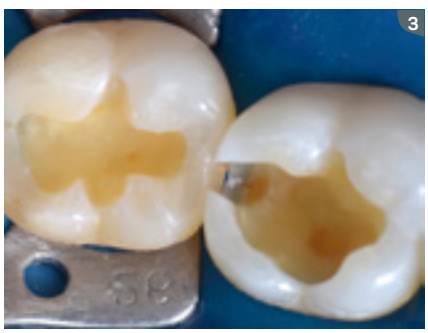
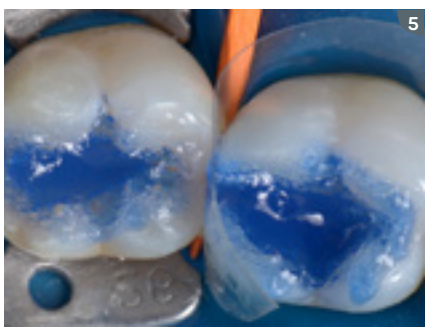
TESTIFIED BY RESEARCHERS APPROVED BY DENTISTS.

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
Ethanol-based solvent.
Solvente a base de etanol.

*MDP performs the chemical bond which, added to the mechanical adhesion already present in the product, gives Ambar Universal a superior adhesion potential.

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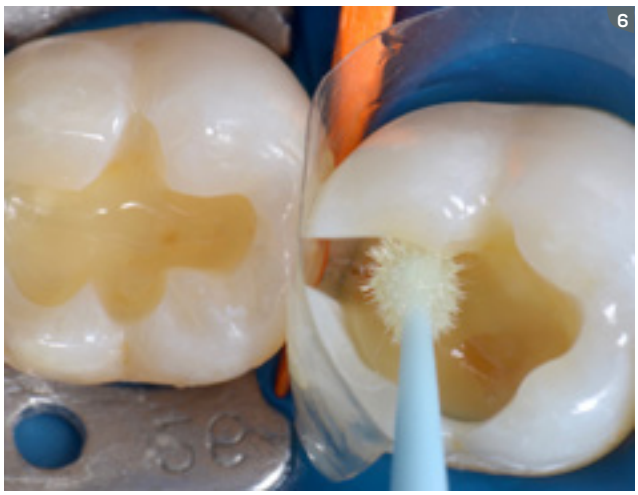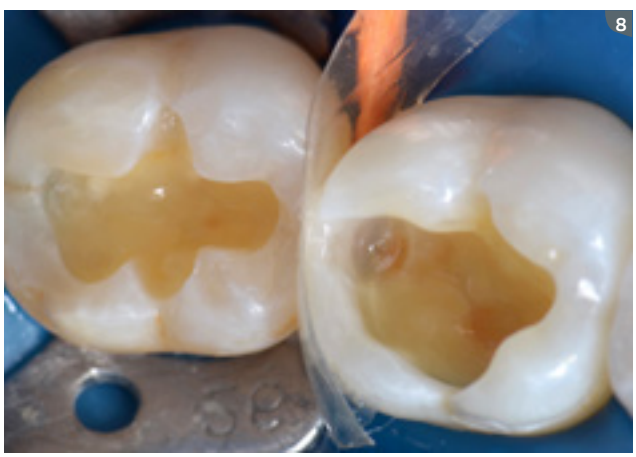
8. Reconstruction of the proximal wall.
8. Reconstrucción de la pared proximal.

9. Restored proximal wall.
9. Pared proximal restaurada.

"Opus Bulk Fill composite meets the needs of any case that performs direct restorations with composites, allowing easy-to-apply, quick and safe protocols. It is the modernity reaching the posterior areas."

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LIGHT-CURING BULK-FILL COMPOSITE
Opus Bulk Fill

RESINA FOTOCURABLE TIPO BULK FILL
Opus Bulk Fill

5th GENERATION LIGHT-CURING ADHESIVE
Ambar

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