The result obtained with whitening with trays was quite satisfactory and did not present any side effect.

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
INITIAL | INICIAL

CASE HIGHLIGHTS | DESTAQUES DEL CASO

It is natural for canines to present greater color saturation than the other teeth. In this case, the canines showed shade A3.

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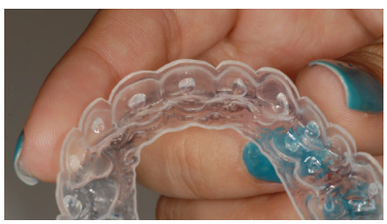It is important to show the patient the appropriate quantity of gel in the tray (one drop per tooth) because some believe that the quantity is related to the quality, which is not true.

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Trimming the tray around the gum edges prevents excess gel from unnecessarily coming into prolonged contact with the gums.

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At the end of the treatment, the comparison with the initial shade guide is beneficial in showing the patient the level of whitening obtained.

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Yellowish teeth.

INITIAL APPEARANCE

Shade A3 for canines and A2 for the other teeth, according to the VITAPAN® Classical shade guide.

CLINICAL/X-RAY EXAM

Good oral health, healthy teeth and gums. It was concluded that the teeth were vital and naturally darkened.

TREATMENT

Supervised take-home whitening with carbamide peroxide at 16% was performed for 28 days, with weekly reviews.

> "The technique of supervised take-home whitening produced excellent results with comfort to the patient, without any side effects."
>
> *DR. ANTONIO SAKAMOTO JUNIOR*

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FGM MATERIALS USED

TAKE-HOME DENTAL BLEACHING GEL (CARBAMIDE PEROXIDE)
Whiteness Perfect 16%

PLATES FOR TRAYS
Whiteness Plates for Custom Dental Trays 1mm

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