EN Microabrasion and whitening as a **conservative** alternative

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EN CLINICAL/RADIOGRAPHIC EXAM

Teeth showing hypoplasia on enamel with whitened and brownish stains in all teeth. When investigating the possible factors for the disease, the patient reported not having a family history of those lesions and that, during deciduous teething she did not show any alterations on the enamel. However, chicken pox and rubella were reported before the age of 6. During the anamnesis, it was not possible to identify intoxication by fluoride during childhood. Nevertheless, because of the characteristics of the lesions and because they affected all the teeth, the exam suggested chronic intoxication through fluoride from birth to the beginning of the eruption of permanent teeth.

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TREATMENT

Microabrasion combined with at-home whitening was suggested as a conservative treatment. Initially, the at-home whitening was carried out with White Class 7.5% for one month. During that period, the professional carried out two microabrasion sessions

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with Whiteness RM for each arch, in intervals of one week. In each appointment, the paste was applied 7 times. After microabrasion, one more month of at-home whitening was carried out, concluding the treatment.

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STEP BY STEP | *PASO A PASO*


1a to 1c. Smile of the patient and intrabuccal aspect evidencing stains on inferior teeth.

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
2. First microabrasion session with absolute isolation installed. To increase the protection of the marginal gingiva, the gingival barrier Top dam was appllied.

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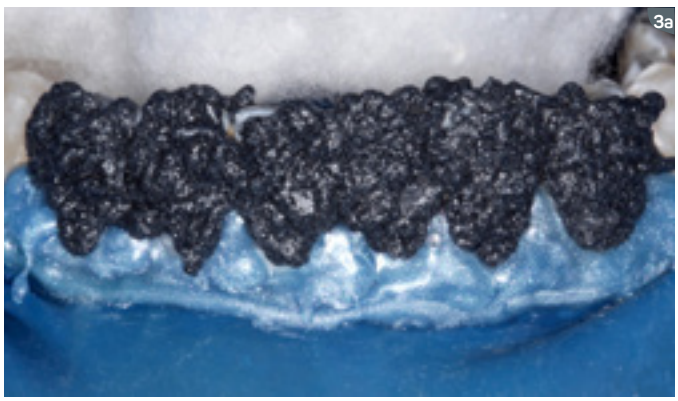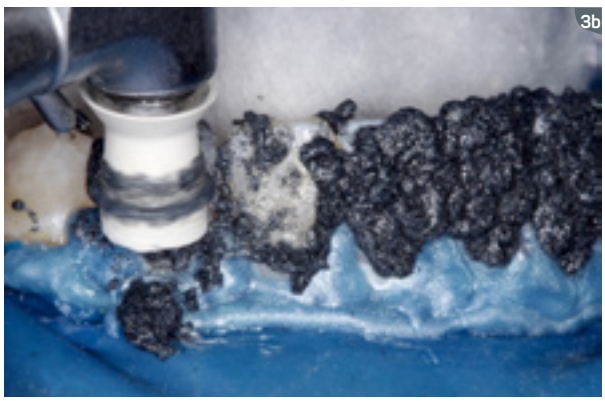
3a and 3b. The paste has to cover the stained surface completely and then the microabrasion can be carried out with rubber cups or plastic spatulas (manually). Each application comprehends the friction for 15 seconds on the stains. In this case, there were 2 sessions with 7 applications each.

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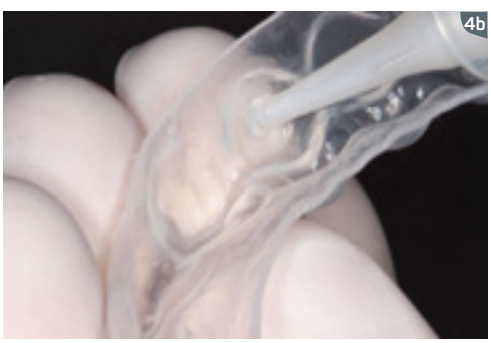4a and 4b. Aspect after the session, the teeth show rugosity and dehydration/demineralization. To accelerate the remineralization process, it is possible to use Desensibilize KF 0.2%, which contains fluoride in its formula. In this case, the gel was applied in a tray; however, it could be done directly at the office with the 2% concentration (Desensibilize KF 2%).

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"I confess that I was impressed with the result obtained after the treatment with microabrasion and at-home whitening. The possibility of removing stains in a simple, practical and efficient way is surprising. All colleagues and professionals I talked to proposed the treatment with ceramic veneers and crowns. Within my own conception of the case, I thought it would be worth a shot to try microabrasion independently of the result, and, to my surprise, it was better than I expected."

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Prof. Dr. Jonas Andrade

5. Aspect of the re-hydrated and remineralized teeth one week after the conclusion of the whitening treatment.

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FGM MATERIALS USED |
MATERIALES FGM UTILIZADOS

STAIN REMOVER BY MICROABRASION
Whiteness RM

TAKE-HOME WHITENING GEL
(HYDROGEN PEROXIDE)
White Class 7,5%

LIGHT-CURING GINGIVAL BARRIER
Top dam

DESENSITIZING GEL
Desensibilize KF 0,2%

PLATES FOR TRAYS
Whitenes Plates for Custom Dental Trays

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